

PRODUCT FEATURES:

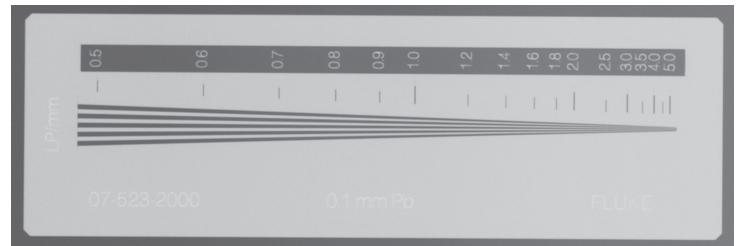
- Perform quick and easy X-ray machine focal spot resolution testing
- Evaluate magnification techniques
- Choose between different versions with different thicknesses and resolution ranges from 1-20 LP/mm

Test Patterns


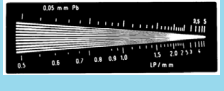
Test patterns are used to determine the resolution of the X-ray machine.

Different number of line pairs per mm (LP/mm) and different thickness of the lead shield is used depending on application and modality.

Our high-precision X-ray test patterns vary from 1 - 20 LP/mm.



TECHNICAL DATA

Model 07-501-2000	1.0-4.8 LP/mm, 16 groups, 0.1 mm Pb, 110x40 mm	
Model 07-523-2000	0.5-5.0 LP/mm, 1 group, 0.1 mm Pb, 157x50 mm	
Model 07-539	1.5-20.0 LP/mm, 1 group, 0.025 mm Pb, 60x30 mm	
Model 07-535	0.6-5.0 LP/mm, 20 groups, 0.05 mm Pb, 50x50 mm	